



Sample Batch Granulator laboratory or small plant use...

The **Granu-Grinder** is a sample batch granulator for laboratory or small plant use. It can be used in the laboratory for granulating small experimental compound batches or used to grind scrap molded parts.

Commonly used for grinding thermo-plastic materials the **Granu-Grinder** is extremely useful for grinding

- elastomers
- composites
- ceramics
- and other various compounds...

The design of the rotor configuration and the low rotational speed of the rotor ensure excellent quality granulate, with consistent size without dust or fines.




It features automatic safety shutoff contact switches for the main housing and recovery trough to ensure safe operations.

- Equipped with
- 1.8kW motor
 - Stationary Cutters
 - 150mm Rotary Cutters

Various screen sizes (Standard 5mm) are available to provide a vast range of granular proportions. Feed opening 5" (127mm) x 9.75" (247mm).

Mounted on a mobile floor stand with four heavy-duty casters.

Granu-Grinder

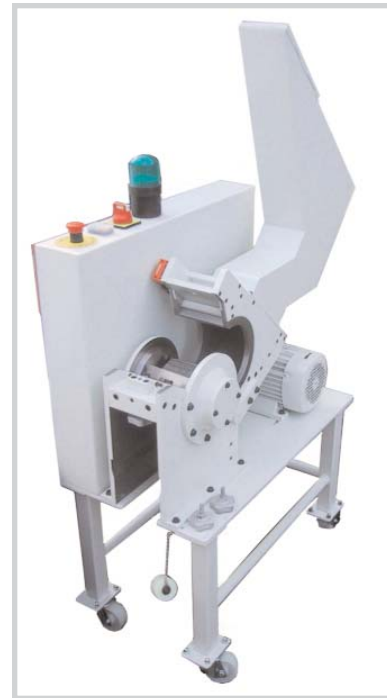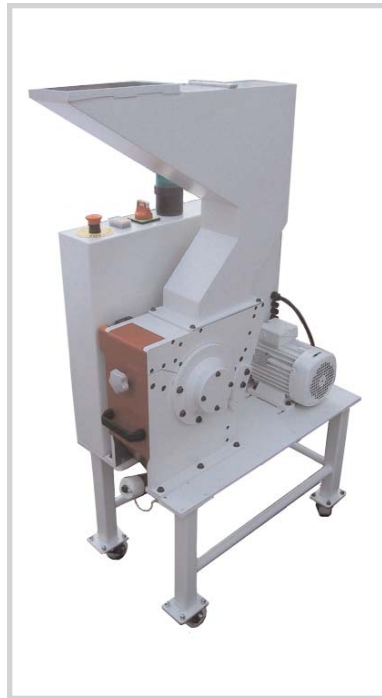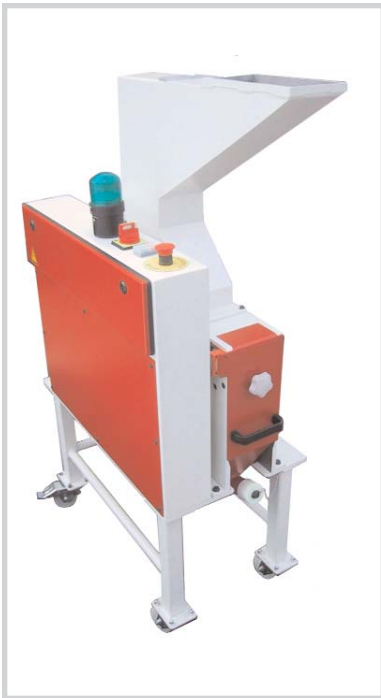

Features and Specifications

C.W. Brabender® • Granu-Grinder

Customer Service

C.W. Brabender® Instruments, Inc. provides unparalleled service and technical support for our customers by employing highly skilled tradesmen, service technicians, and an experienced sales force.

We have a modern application laboratory located at our National Headquarters to benefit the interests of our customers. An experienced technician shall attend to the specific needs of each and every individual, and shall remain present throughout the entirety of the test and trial periods in order to assist in the customer's quest for desired results.

To arrange for a personal demonstration of the Granu-Grinder, contact the technical staff at C.W. Brabender® .

Discuss what WE can do
for YOU...



Technical Data

Power 1.8kW

Motor Speed 1720 RPM

Feed Opening 127mm x 247.7mm • 5" x 9.75" (W x L)

Power Requirements 230 VAC, 60 Hz, 7.8 Amp, 3-Phase

equipped with (2) stationary cutters
(3) 150mm rotary cutters
(1) 5mm screen (hardened)
built in safety contact switch



C.W. Brabender®
INSTRUMENTS, INC.

C.W. Brabender® Instruments, Inc. • 50 East Wesley Street, South Hackensack, NJ 07606
Tel: 201.343.8425 • Fax: 201.343.0608
E-Mail: cwbi@cwbrabender.com • www.cwbrabender.com

